Dutchblue.com Nederland

Markt 30a
2611 GV Delft
Nederland

T 00 800 2611 2611 (gratis) of 015 3030 025
F 084 883 66 06
E info@dutchblue.com
W www.dutchblue.com

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Square de Meeûs 37
1000 Brussel
België

T 00 800 2611 2611 (gratis) of (0032) 026 699 663
F (0032) 02 403 05 55
E info@dutchblue.com
W www.dutchblue.com

Dutchblue.com Deutschland

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Dutchblue.com United Kingdom

Park House
111 Uxbridge Road
Ealing, London, W5 5TL

T 00 800 2611 2611 (free) or 02 033 977 242
F 02 033 977 524
E info@dutchblue.com
W www.dutchblue.com

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
REF 5000
IFU 5000 – 13M2

ESCAPIES evacuation mat

The ESCAPIES evacuation mat is placed between the mattress and the bed slats, allowing nursing home residents to be evacuated by a single staff member. Bed-ridden residents lying on an ESCAPIES evacuation mat can be slid across a variety of surfaces and through f.e. narrow staircases.

Instructions for use

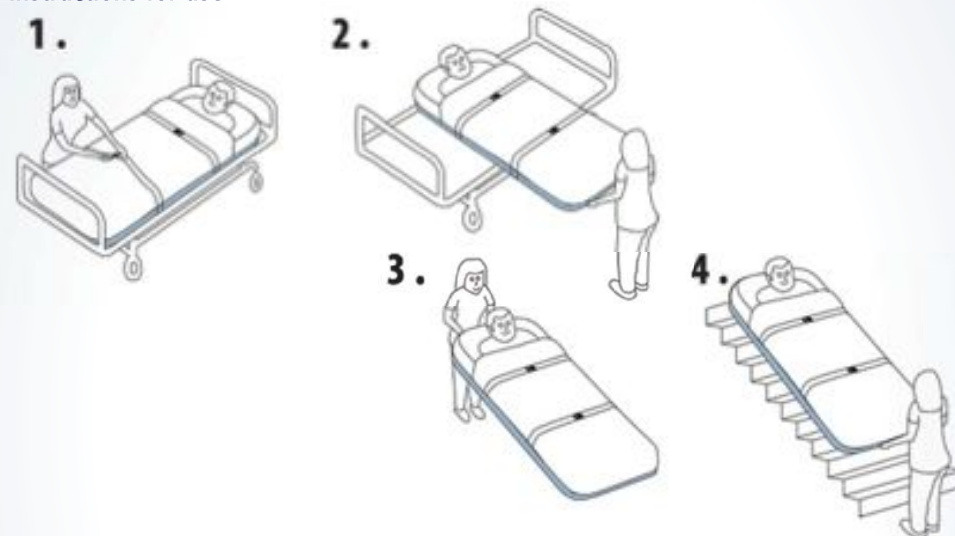- Remove the packaging. Do not use scissors, knives or other sharp implements when opening the packaging.
- Insert the evacuation mat between the mattress and the bed slats of the nursing bed as a standard.
- Attach the mat to the nursing bed mattress using the elastic bands at the corners.
- In the event of an emergency securely fasten the nursing home resident to the evacuation mattress using the belts.
- Use the grips at the head and/or foot of the evacuation mattress to slide the mattress (with the resident on it) across the floor into safety.

Instructions for cleaning

The evacuation mat can be cleaned with a soft cloth and mild detergents or clinical cleaning products. Never use disinfection agents, phenol or solvents.

PLEASE NOTE:

- Do not allow to come into contact with open fire.
- Avoid contact with sharp objects.

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