Dutchblue.com Nederland

Markt 30a
2611 GV Delft
Nederland

T 00 800 2611 2611 (gratis) of 015 3030 025
F 084 883 66 06
E info@dutchblue.com
W www.dutchblue.com

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Square de Meeûs 37
1000 Brussel
België

T 00 800 2611 2611 (gratis) of (0032) 026 699 663
F (0032) 02 403 05 55
E info@dutchblue.com
W www.dutchblue.com

Dutchblue.com Deutschland

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Dutchblue.com United Kingdom

Park House
111 Uxbridge Road
Ealing, London, W5 5TL

T 00 800 2611 2611 (free) or 02 033 977 242
F 02 033 977 524
E info@dutchblue.com
W www.dutchblue.com

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User Manual – Viscoelastic mattress

REF 3020 – 3030
IFU 3020 3030 – 13M2

dutchblue.com

DREAMIES viscoelastic mattress

DREAMIES nursing mattresses provide a high degree of comfort. They mould themselves to the user's contours while offering optimum support of the entire body.

The top layer of the mattress is made from a special type of viscoelastic memory foam that facilitates an even distribution of pressure and body weight. The top layer of PREMIUM mattresses is made from a sturdier layer of viscoelastic memory foam that facilitates an even better distribution of pressure and body weight. This improves circulation of the blood and the supply of oxygen to the tissues. The sturdy bottom layer ensures that the entire body is supported to the greatest extent possible. Effective for decubitus prevention.

Composition

Mattress


- 4 cm viscoelastic memory foam
50 kg/m³ (Basic) 85 kg/m³ (Premium)
- 8 cm polyurethane foam

Mattress cover

Top	100% PU
Bottom	100% polyester

Basic REF 3020

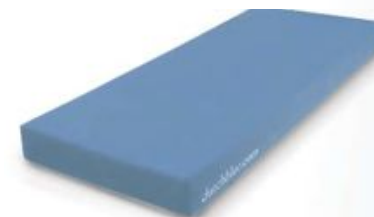Premium REF 3030

Instructions for use

- Remove the packaging. Do not use scissors, knives or other sharp implements when opening the packaging.
- Unroll the mattress and leave unused for one day, in order for the memory foam to fully expand.

Instructions for cleaning

The mattress cannot be washed, but can be aired out. The cover can be cleaned with a soft cloth and mild detergent or clinical cleaning products. Never use disinfection agents, phenol or solvents.

PLEASE NOTE:

- The typical odour of visco-elastic foam will disappear in time.
- Do not allow to come into contact with open fire.
- Avoid contact with sharp objects.

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